



**Cricath[®], the smart choice
for cricothyrotomy**







Cricath[®], the smart choice:

- › small lumen – 2 mm ID
- › of flexible, high grade Teflon[®]
- › and an innovative tilting mechanism
- › to reduce the risk of kinking

Cricath[®] is a flexible needle cricothyrotomy catheter that enables a quick access to the trachea. This application asks for the best materials. The lumen of Cricath[®] is made of a special, high grade type of Teflon[®] (PFA). The patented tilting mechanism enables a proper fit of the lumen into the trachea, also in case of an abnormal neck anatomy. Together with the elevated fixation point of the lumen and the enlarged, smoothly curved opening in the flange around the lumen, this provides for additional flexibility, reducing the risk of kinking. With the standard catheter over needle procedure for cricothyrotomy, Cricath[®] can be placed by insertion through the cricoid membrane.

cricath[®], the smart choice

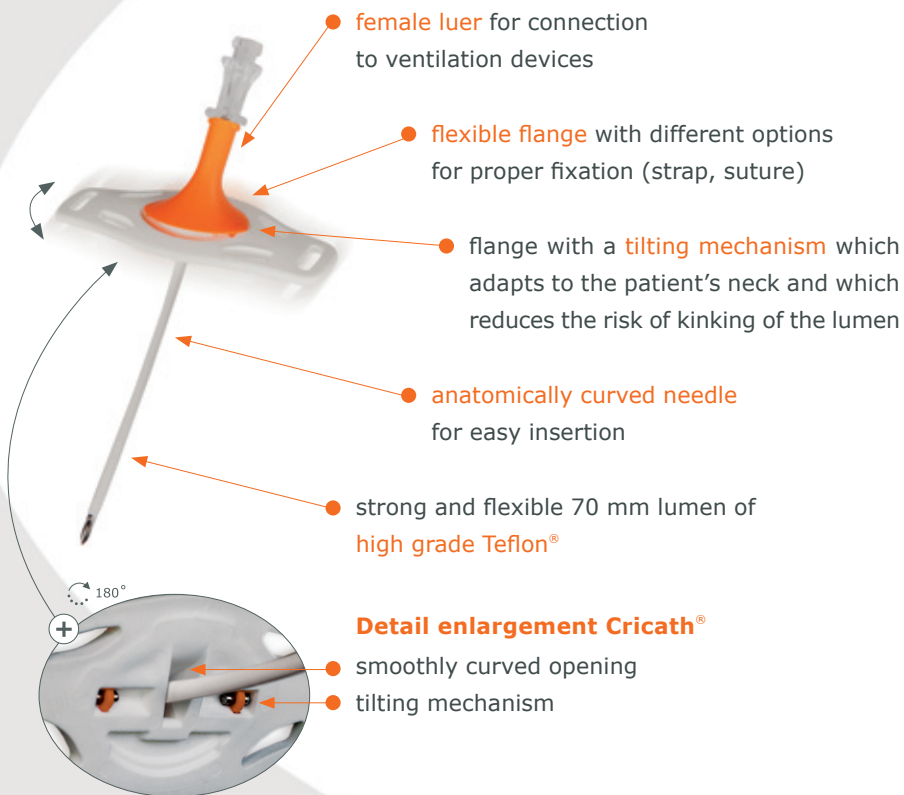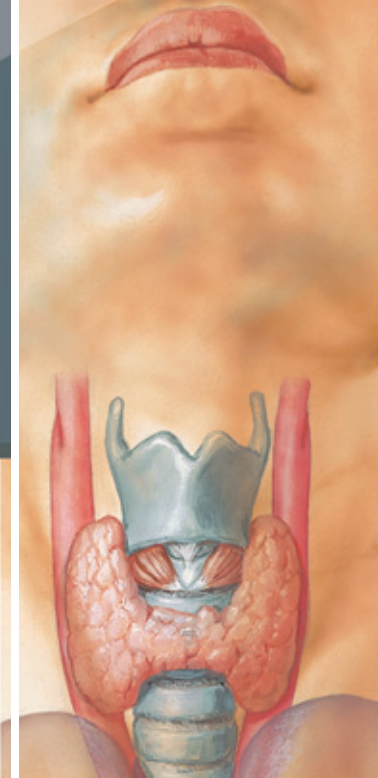

Cricath® is protected by a patent application and a design registration.

»
**quick
access to
the trachea**

ce for cricothyrotomy



*Paola Vasters
product manager*

*Dirk van Asseldonk
CEO*

*Caroline Abeln
customer services*



Development

Cricath[®] is developed in the Netherlands by Dolphys Medical, the same company that developed Ventrain[®]. For more information about Cricath[®] and our other products please visit our website or contact us directly.

Ordering information Cricath[®]

- packed with syringe and neck strap
- single use
- sterile
- art. code 2618400

Website

www.dolphys.com/cricath

Your distributor

Dolphys Medical B.V.

Info Paola Vasters
Caroline Abeln

P De Lismortel 31
5612 AR Eindhoven
The Netherlands

T +31 (0)40 751 60 20
E cricath@dolphys.com
I www.dolphys.com